RYZRAIRBOOTS

User Manual

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NOW YOU
ARE A
RYZR

SAFETY AND WARNING INFORMATION

The RYZR Airboots are a modern massage device with electrical components. They have been designed to operate safely and reliably. However, please note that improper use or maintenance may lead to injury or technical issues.

To ensure safe use of your RYZR Airboots, please read the following safety information carefully and in full before first use.

Contraindications

The RYZR Airboots are intended for use by healthy adults and are particularly suitable for recovery after physical activity or a long day. If you are currently undergoing medical treatment or are unsure whether the device is suitable for you, please consult your doctor before use – especially if you have recently had surgery or have any pre-existing medical conditions.

Do not use the device in the following cases:

- Respiratory conditions such as pulmonary oedema or chronic pulmonary insufficiency
- Heart failure, severe cardiac conditions or cardiac pacemakers
- Thrombosis, blood clotting disorders, inflammatory venous conditions
- Wounds, infections or inflammatory skin conditions in the application area
- Arterial circulatory disorders or advanced arteriosclerosis
- Pregnancy – no clinical studies on use during pregnancy are available
- Recent skin grafts
- Kidney disorders accompanied by oedema
- Implants or artificial joints in the direct application area
- Tumours, fractures, severe pain or numbness in the legs

General Notes on Use

- Use the RYZR Airboots exclusively on the legs.
- Keep the RYZR Airboots out of reach of children. Please ensure that children do not use or play with the RYZR Airboots unsupervised.
- Persons with limited physical or cognitive abilities should only use the RYZR Airboots under supervision.
- Do not consume alcohol before use.
- In case of symptoms such as pain, shortness of breath, dizziness or bruising: stop use immediately and seek medical advice.
- Adopt a comfortable, relaxed position – slightly elevating the legs may enhance the effect.
- We recommend starting at the lowest pressure level and increasing gradually.

RYZR Airboots comply with applicable EU safety, health, and environmental protection requirements.

Additional safety instructions

- Do not make any technical modifications to the device. Modifications may cause safety hazards or malfunctions and may void the warranty.
- If discomfort occurs during use, you can activate the emergency relief valve at any time. This immediately releases the pressure, allowing you to safely remove the Airboots. Afterwards, set the valve back to “OFF”.
- Do not use the RYZR Airboots near flammable gases, heat sources or strong electromagnetic fields.
- Protect the RYZR Airboots from moisture, corrosive liquids and dust. No liquid must enter the device.
- Do not place sharp objects (e. g. needles) near the RYZR Airboots.
- Keep the RYZR Airboots away from children and pets – risk of suffocation and strangulation!
- Regularly remove dust and lint from the controller and hose connections – preferably with a dry cloth.
- Do not use the RYZR Airboots while sleeping.
- Check the condition of your RYZR Airboots before each use. Do not use them if you detect leaks or visible damage.
- Ensure that the zips are closed and any accessories on feet or ankles are removed before starting the RYZR Airboots.
- Do not use the RYZR Airboots together with other electronic devices or medical equipment. This may cause malfunctions.
- In the event of technical problems or malfunctions, switch off the RYZR Airboots and contact RYZR Customer Service.
- Do not drop the RYZR Airboots or subject them to impacts.
- Ensure suitable environmental conditions during use (5–40 °C, max. 90 % humidity, max. 3,000 m altitude).
- Do not machine wash the RYZR Airboots. If necessary, gently wipe the surface with a damp cloth and mild detergent. Avoid getting liquid into the air connections.
- The RYZR Airboots must not be wrung out, ironed, tumble-dried or dried with other heat sources. Do not use bleach.

Safety Instructions for Lithium-Ion Batteries

- Use only the batteries supplied with your RYZR Airboots.
- Do not open, disassemble or deform the battery – risk of short circuit, fire or explosion.
- Keep the battery away from fire, heat sources, direct sunlight and water.
- Do not store batteries together with metal objects (e. g. keys, chains, hairpins).
- If the battery swells, overheats, leaks, becomes deformed or emits an unusual odour, stop using it immediately and contact RYZR Customer Service.
- If battery fluid leaks, avoid skin and eye contact. If the fluid contacts skin, rinse thoroughly with water immediately and seek medical advice if symptoms persist. If it contacts eyes, do not rub, rinse immediately and thoroughly, and seek medical attention urgently – the fluid may cause severe eye injuries or blindness. Remove any clothing that has come into contact with the fluid immediately and rinse the skin thoroughly.
- If you do not use your RYZR Airboots for an extended period, recharge the battery approximately every three months to preserve its lifespan.
- Always keep batteries out of reach of children. Swallowing or improper handling may cause serious injury or be life-threatening.

PRODUCT

PRODUCT SPECIFICATIONS

Product Name	RYZR Airboots (Compression Massage Device)
Model	S: A20027714 M: A20027715 L: A20027716 XL: A20027717
Pressure Range	Adjustable between 40–180 mmHg
Rated Power	25 W
Rated Voltage	7,3 V 18,98 Wh
Input Voltage	5 V 2 A
Battery Capacity	2600 mAh
Massage Programmes	M1 Fatigue Relief M2 Training Activation M3 Workout Recovery M4 Overall Relaxation M5 Deep Care
Time Settings	15–60 min in 5-minute increments
Pressure Levels	30 levels, 5 mmHg per level
Battery/Accumulator	Not permanently installed

SIZE CHART AIRBOOTS

	L × W × H (cm)
S	84.0 × 38.0 × 14.0
M	94.0 × 38.0 × 14.0
L	102.0 × 38.0 × 15.0
XL	112.0 × 40.0 × 15.0

Control unit with display

Emergency regulator

High-quality zipper

Replaceable patch

SCOPE OF DELIVERY

1 pair of RYZR Airboots including storage bag
 2 charging cables
 2 batteries including storage pouch

Documents included in delivery:
 1 Quick Start Guide
 1 Customer Service contact card

Rechargeable batteries / 2600mAh


Charging cable

Storage bag

MAIN STRUCTURE OF THE DEVICE


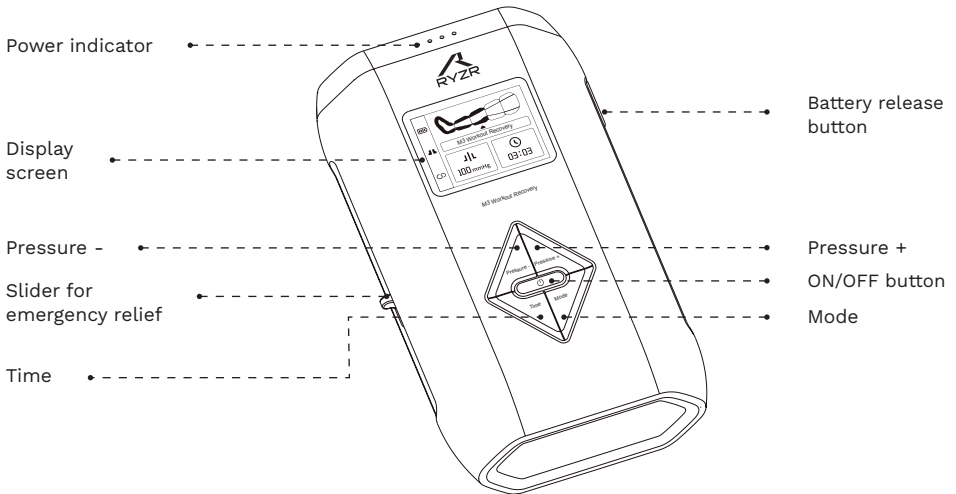
COMPONENTS AND FUNCTIONS

App connection

Chamber status Area intensification

Battery status


Coupling

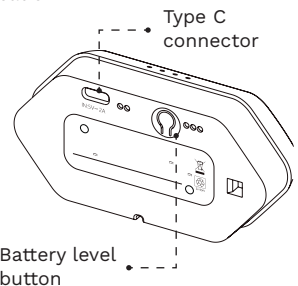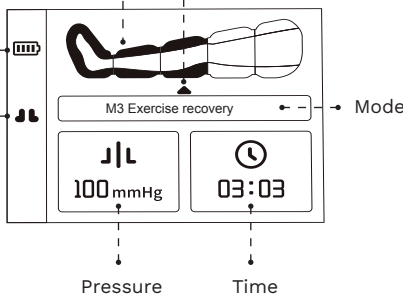OPERATION


Inserting the Battery

Carefully insert the charged batteries into the controller. Ensure they are fully engaged – you will hear a distinct “click” when correctly seated.

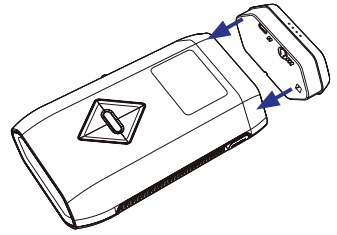Removing the Battery

First end the session and switch off the device. Then press the battery release button marked “OPEN” on the side of the controller. The battery will spring out slightly and can be removed for charging or stored in the supplied pouch.

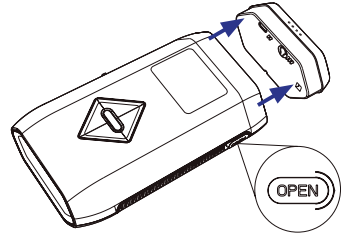Before First Use


Before using your RYZR Airboots for the first time, prepare them carefully to ensure everything works smoothly:

- Fully charge the batteries: Charge both batteries completely – ideally for approximately three hours. Use the supplied USB-C cable and connect it to a power bank or a mains adapter with 5 V / 2 A output. If your adapter delivers less than 2 A, charging time will be extended accordingly.
- Check the charge status: During charging, the LEDs on the battery flash. Once all four LEDs are lit continuously, the battery is fully charged. If only the first LED is lit, the battery level is very low – please charge promptly.
- Insert the batteries: Carefully insert the charged batteries into the controller.
- Inspect: Check all components shortly before use. Look for damage, twists or leaks. If anything appears unusual, do not use the RYZR Airboots and contact RYZR Customer Service (support@ryzr.de).

Application

- Do not use the RYZR Airboots near flammable gases, heat sources or strong electromagnetic fields.
- Put on the Airboots: Pull the RYZR Airboots over both legs and ensure they fit correctly. Make sure the zips are fully closed.
- Get comfortable: Adopt a relaxed lying or sitting position. Slightly elevating the legs can additionally support lymph flow and circulation.
- Switch on the device: Briefly press the power button. The controller display will light up – the RYZR Airboots are now in standby mode.
- Select a mode: Use the MODE button to choose one of the five preset programmes. They differ in sequence and intensity.
- Set the pressure: Adjust the pressure levels using the “PRESSURE +” and “PRESSURE –” buttons. You can also change the pressure during an active session. The RYZR Airboots will pause briefly and then resume the massage at the new pressure level.
- Set the duration: Use the “TIME” button to set the application duration in 5-minute increments – between 15 and 60 minutes.
- Start the massage: Press the power button again to start the session. If you wish to end the massage early, press the same button again.
- Use the emergency function: If you experience discomfort during use, slide the “Emergency Pressure Relief” switch to “ON”. The air pressure is released immediately and you can safely remove the Airboots. Set the switch back to “OFF” when you wish to use the device again.

- **Switching off & deflate function:** Once the set time has elapsed, the RYZR Airboots automatically return to standby mode. If no input is made for more than 60 seconds, the RYZR Airboots switch off automatically. You can manually switch off at any time by holding the power button for two seconds. The air in the RYZR Airboots is released automatically (approx. 2 minutes). To skip the deflation process, hold the power button again for two seconds to switch off immediately.

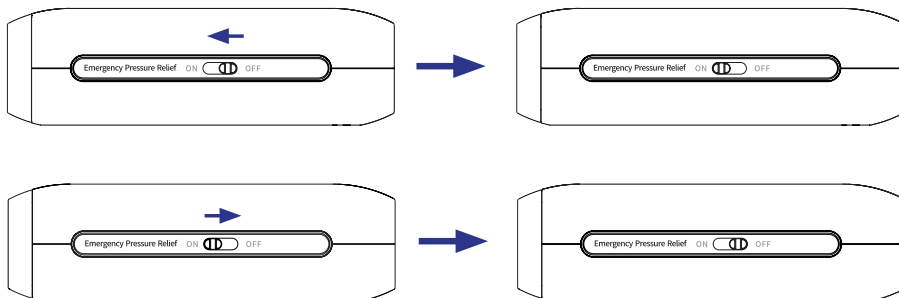INDIVIDUALLY DEACTIVATING MESSAGE ZONES

Your RYZR Airboots feature five individually activatable massage chambers, ranging from the foot to the upper thigh. By default – regardless of the selected mode – all five chambers are activated. If required – e.g. on medical advice or for targeted relief of specific areas – individual chambers can be deactivated before starting a session.

In standby mode, hold the “PRESSURE +” button for approximately one second. The controller switches to the graphical chamber selection view. Use the “PRESSURE +” and “PRESSURE –” buttons to select the desired active chambers.

Activated chambers are shown as light grey symbols on the display. Deactivated chambers appear correspondingly darker. To confirm your selection, briefly press the start button – you will return to the home screen.

Note: The last selected configuration is saved and remains stored even after the RYZR Airboots are switched off. Changes during an active session are not possible. End the current session to make a new selection.

EXTENDING MESSAGE TIME FOR SELECTED ZONES

In “Workout Recovery” mode, you can target specific massage zones for more intensive treatment. The selected zone is massaged twice additionally within one session – at the same pressure level. To do this, hold the “PRESSURE –” button for approximately one second. The display shows a graphical selection view where you can choose the desired zone using the “PRESSURE +” and “PRESSURE –” buttons. The selected zone is highlighted in the view. Confirm your selection by briefly pressing the start button to return to the main screen. This allows you to give specific muscle groups extra attention without changing the overall sequence.

This function is only available in “Workout Recovery” mode and can only be set before starting a session. Changes during an active massage are not possible.

Figure 1

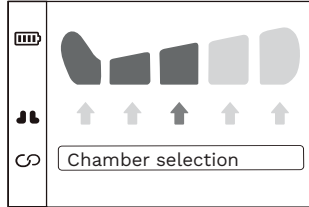
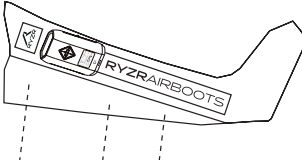
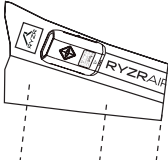
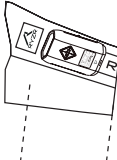
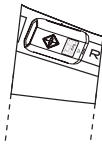
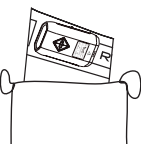
Figure 2


STORAGE

After use, you can fold your RYZR Airboots compactly and store them in the supplied storage bags.

<p>Step 1: After decompression, fold flat.</p>	<p>Step 2: Fold the two ends together.</p>	<p>Step 3: Fold in half again.</p>
		
<p>Step 4: Fold the ends completely together.</p>	<p>Step 5: Open the storage bag and place the RYZR Airboots inside.</p>	<p>Step 6: Finally, tie the storage bag closed and set aside.</p>
		

MODE DESCRIPTIONS

Fatigue Relief and Training Activation Mode:

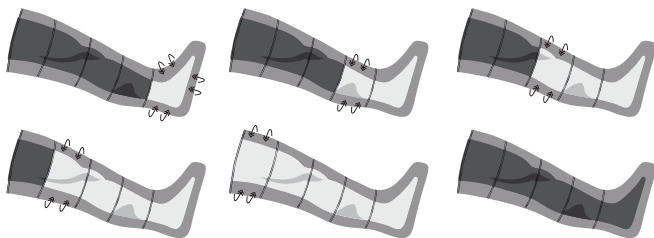
Fatigue Relief:
(Relieves and revitalises heavy legs)

In “Fatigue Relief” mode, the massage begins at the foot and moves upward towards the upper thigh. The compression chambers are activated sequentially until all zones are evenly pressurised. A complete pressure release then follows. The aim is to promote venous return, reduce fluid accumulation in the legs and relieve the general feeling of heaviness.

Training Activation
(Prepares the muscles for exertion)

This mode is designed to prepare the body before physical activity. The massage follows the same ascending sequence as “Fatigue Relief” – from bottom to top. Compared to “Fatigue Relief” mode, the pressure is held for longer, specifically preparing the muscles for upcoming exertion.

Workout Recovery Mode:

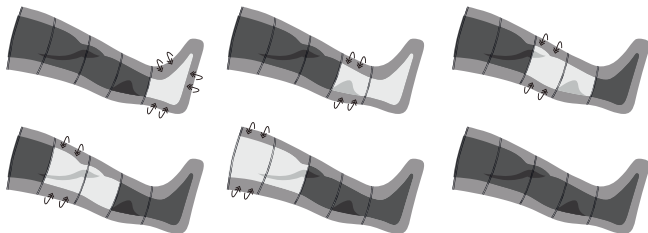Workout Recovery:
(Supports recovery after physical activity)


This mode is specifically designed for recovery after physical activity and offers a deep-acting, circulation-promoting massage for sustained relief of muscle soreness and muscular fatigue.

The sequence follows a smooth, overlapping rhythm: first, the lowest chamber in the foot area is filled with air and the pressure is briefly held. The next chamber then begins to fill while the previous one releases. This alternating sequence continues until the uppermost chamber is reached.

Throughout the session, two chambers are always active: one holds pressure while the next is already being filled. This creates a flowing, gliding pressure pattern along the leg – comfortable, effective and soothing.

In this mode, you can target specific zones for more intensive treatment or exclude them entirely (see pages 9 and 10).

Overall Relaxation Mode:

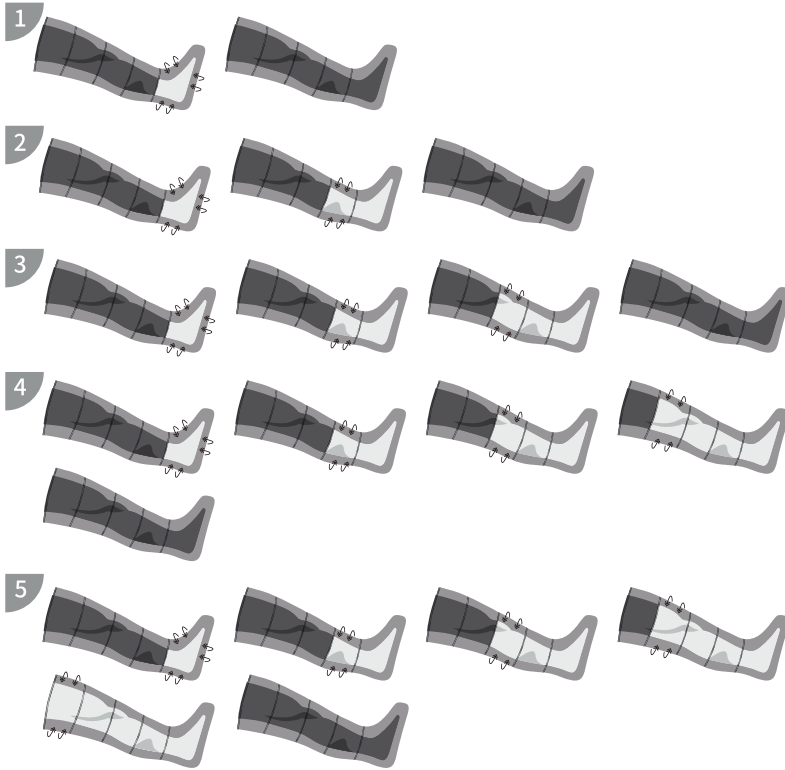Overall Relaxation:
(Gentle massage for
all-round soothing
relaxation)


The “Overall Relaxation” mode uses a pleasant, progressive massage sequence. Initially, only the lowest chamber in the foot area is activated. It fills with air, the pressure is briefly held and then fully released. This process is repeated once before additional chambers are added.

In the next cycle, the first chamber is filled and held again – then the second chamber also begins to fill. After a brief pressure phase, both chambers are released simultaneously.

With each further repetition, the sequence extends by one additional chamber – first the third, then the fourth and finally the fifth. The pressure is briefly held before all activated chambers are released simultaneously. Once all five chambers have been cycled through, the sequence begins again.

This calm, progressive rhythm creates an even massage along the entire leg. The application supports circulation, relieves tissue pressure and contributes noticeably to recovery – ideal for relaxed evenings or recovery after long periods of standing or sitting.

Deep Care Mode:

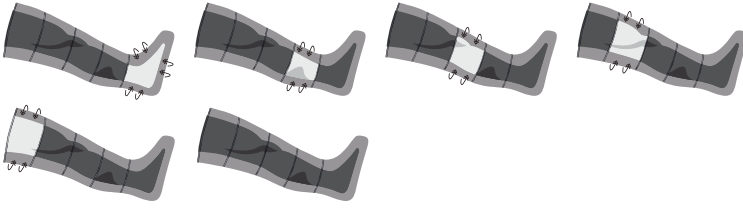Deep Care:
(Targets deep-seated
tension)

The “Deep Care” mode offers a particularly targeted massage to release deep-seated tension. Each chamber is activated individually to apply focused pressure to the respective muscle group. The pressure is briefly held and gently released before the next chamber follows.

TROUBLESHOOTING

Problem	Cause	Solution
The RYZR Airboots will not switch on	The battery is completely discharged	Fully charge the batteries before using the RYZR Airboots again
The batteries will not charge	The charging cable is not connected correctly	Check the connection between the controller and the charging cable
	The charging cable is damaged	Use a functioning USB-C cable
	The cable is incompatible	Use a mains adapter or power bank with a USB-A port and sufficient output
The RYZR Airboots switch off during use	The battery is empty / The air connection is interrupted	Fully charge the batteries / Check that all hoses and connections are properly connected
 The air chambers do not fill correctly	Air circulation is disrupted	Ensure that all air connections are secure and not obstructed
Error message or unusual behaviour	Technical fault, e.g. pressure sensor defective	Please contact Customer Service: support@ryzr.de

For any other issues, please contact: support@ryzr.de

CUSTOMER SERVICE

If you have any questions or need assistance, please do not hesitate to contact us! Our Customer Service team will be happy to help you quickly and straightforwardly. Simply send us an email at support@ryzr.de

If your product is defective or you wish to send it in for inspection, please contact us immediately and have your order number ready.

RoHS DECLARATION OF CONFORMITY

Our products comply with the Directive on the restriction of the use of certain hazardous substances (2011/65/EU) (RoHS 2.0). The raw materials and components of our products comply with the restrictions for the following toxic and hazardous substances.

Cadmium (Cd): <100 ppm	Lead (Pb): <1,000 ppm
Mercury (Hg): <1,000 ppm	Hexavalent chromium (Cr VI): <1,000 ppm
Polybrominated biphenyls (PBB): <1,000 ppm	Polybrominated diphenyl ethers (PBDE): <1,000 ppm
Bis (2-ethylhexyl) phthalate (DEHP): <1,000 ppm	Benzyl butyl phthalate (BBP): <1,000 ppm
Dibutyl phthalate (DBP): <1,000 ppm	Diisobutyl phthalate (DIBP): <1,000 ppm

WEEE DIRECTIVE

Proper Handling of This Product	
	<p>This symbol indicates that the product must not be disposed of with normal household waste. To prevent potential damage to the environment or human health through uncontrolled waste disposal, it must be recycled responsibly to promote the sustainable reuse of material resources. To return a used device, please use the return and collection system or contact the retailer from whom you purchased the product.</p>

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RYZR AG